

The NorPix Sequence File Format (.seq)

NorPix (C) 2013

Sequence files

Overview

The sequence files created by StreamPix/TroublePix/Hermes while recording (.seq extension) use the NorPix Sequence File Format explained here.

A sequence file is made of a header section located in the first 1024 bytes. The header contains information pertaining to the whole sequence: image size and format, frame rate, number of images etc.

Following the header, each images is stored and aligned on a 8192 bytes boundary. Please note that only the uncompressed sequence format is documented here, compressed sequences are handled in a different way.

Usually, pixels in the images are stored for top-left to bottom-right corner. Immediately following the image data comes 8 bytes, containing the absolute timestamp at which the image has been grabbed. The first 4 bytes hold the date and time, then 2 bytes for the milliseconds and the last 2 bytes are the microseconds.

Example

Here is an example for a sequence of 10 images of 640 x 480 pixels in 8 bit monochrome in which the first image in the sequence file is at an offset of 8192 bytes (no metadata).

Read 640 x 480 or 307200 bytes to get all the image pixels. Then read the next 32 bit (4 bytes) to get the timestamp in seconds, formatted according to the C standard *time_t* data structure (32-bit). Read the next 16 bit (2 bytes) as an unsigned short to get the millisecond precision on the timestamp, then read the last 2 for the microseconds.


The image offset for image at index *i* will be : $8192 + (i * TrueImageSize)$

The timestamp information is located at : $8192 + (i * TrueImageSize) + ImageInfo.ImageSizeBytes$.

More

If you need more information or help retrieving data from our sequence files please contact us at support@norpix.com.

Sequence File Sections

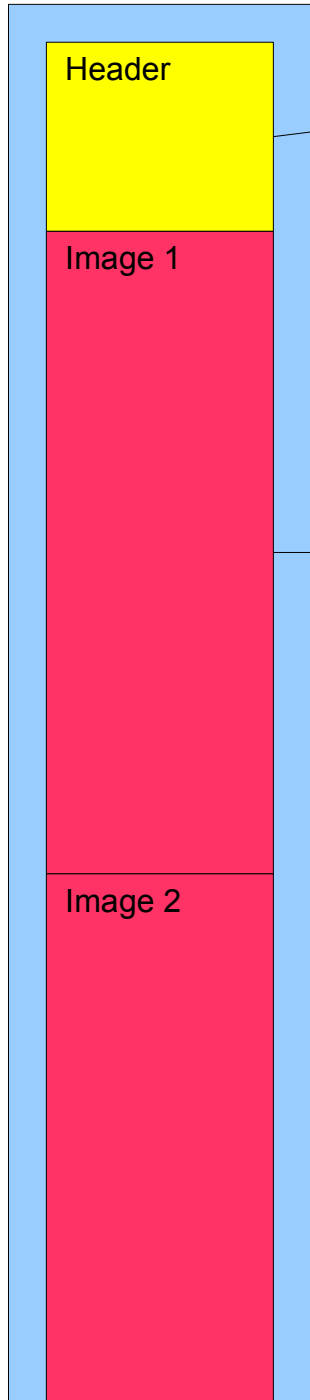

Header


The sequence header hold global information about the sequence file. It has a fixed size of 1024 bytes.

Image 1

In uncompressed sequences, each image section has a fixed size defined in `Header.TrueImageSize`. Each image section holds exactly one image and its time stamp.

Image 2

Sequence File Sections

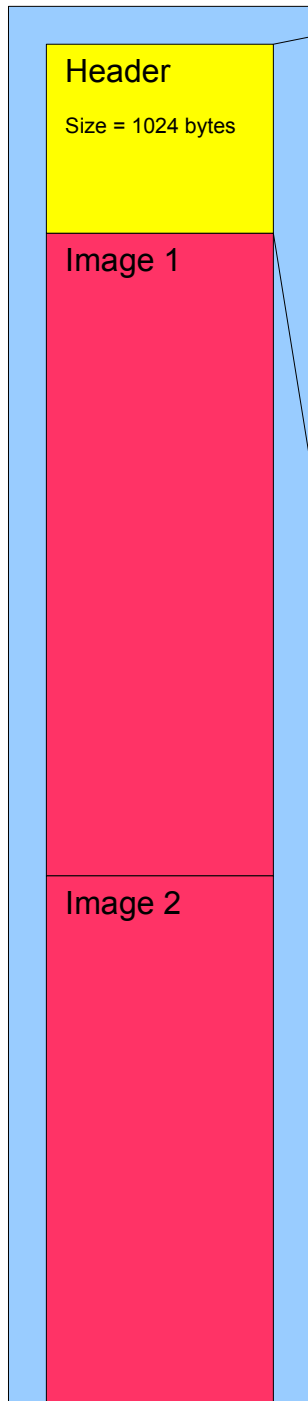

Header


Size = 1024 bytes

Image 1

Image 2

Name	Content	Offset [size]
long MagicNumber	Always 0xFEED	0 [4]
wchar_t Name[12]	Always "Norpix seq\n"	4 [24]
long Version	Sequence Header Version	28[4]
long HeaderSize	Always 1024	32 [4]
BYTE Description[512]	User description	36 [512]
CImageInfo ImageInfo	See belows for a description of the CImageInfo struct	548 [24]
unsigned long AllocatedFrames	Number of frames allocated in the sequence	572 [4]
unsigned long Origin	Should be 0 if not Pre/Post recorded	576 [4]
unsigned long TruelImageSize	Number of bytes between the first pixel of each successive images	580 [4]
double FrameRate	Suggested Frame rate for playback (in fps)	584 [8]
long DescriptionFormat	The content of "Description" 0-unicode 1-ascii 2-data	592 [4]
ULONG ReferenceFrame	ReferenceFrame index (0 if none)	596 [4]
ULONG FixedSize	Fixed size for compressed sequences (0 if none)	600 [4]
ULONG Flags	NorPix reserved flags	604 [4]
long BayerPattern	Bayer pattern used	608 [4]
long TimeOffsetUS	Time offset applied to each image timestamp	612 [4]
long ExtendedHeaderSize	Size of the extended header (if any)	616 [4]
eHCompression CompressionFormat	The compression used	620 [4]
long ReferenceTime	Custom Reference Time (time_t format)	624 [4]
ushort ReferenceTimeMS	Custom Reference Time (milliseconds part)	628 [2]
ushort ReferenceTimeUS	Custom Reference Time (microseconds part)	630 [2]
ulong H264GOP	Group of Picture value if H264 compression is used	632 [4]
ulong H264Bitrate	Bitrate if H264 compression is used	636 [4]
BYTE Padding[384]	Unused bytes, reserved for future uses.	648 [376]

Sequence File Sections

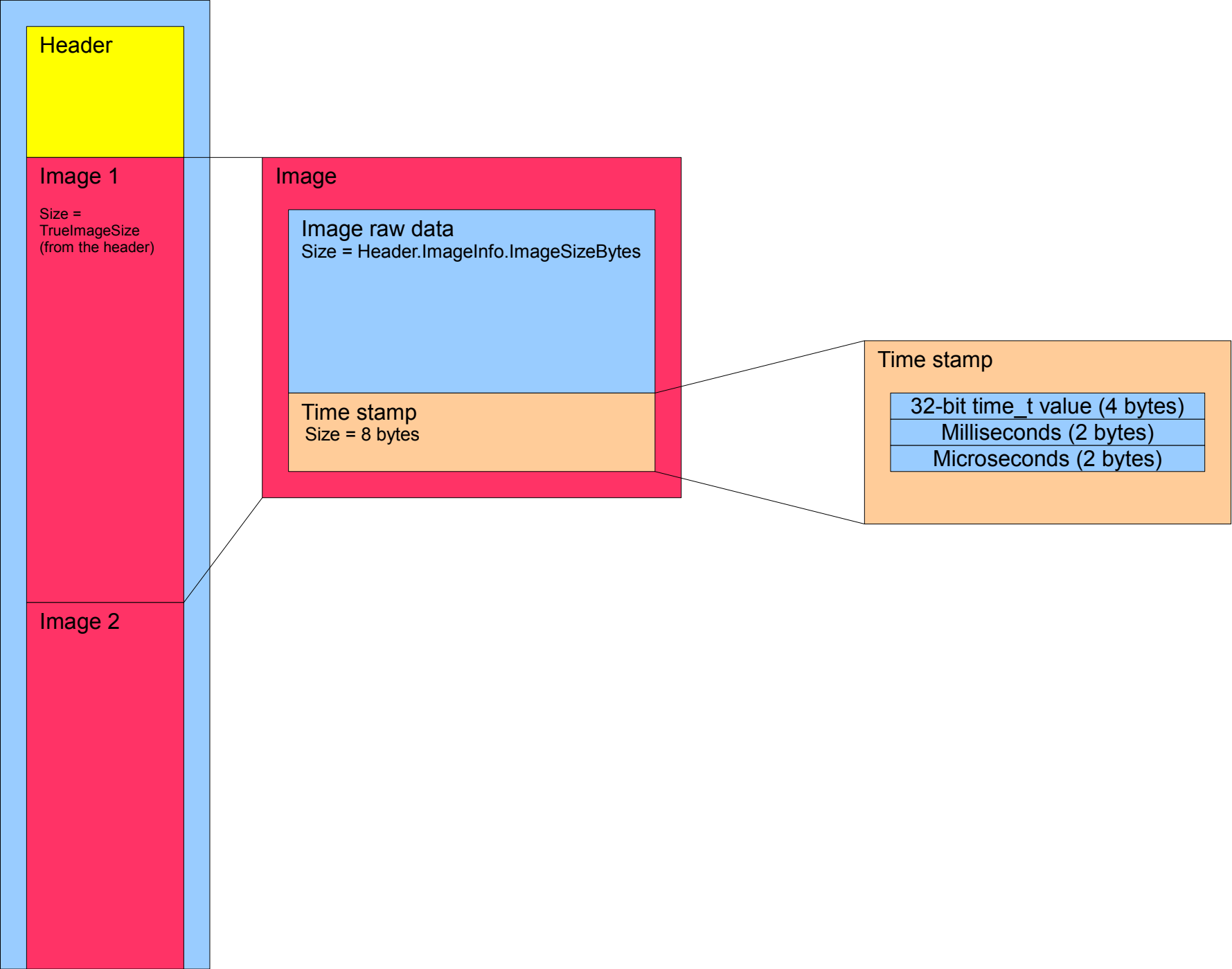

Additionnal Structures and Enums (1/2)

```
struct CImageInfo
```

```
{
    unsigned long ImageWidth;           //Image width in pixel
    unsigned long ImageHeight;          //Image height in pixel
    unsigned long ImageBitDepth;        //Image depth in bits (8,16,24,32)
    unsigned long ImageBitDepthReal;    //Precise Image depth (x bits)
    unsigned long ImageSizeBytes;       //Size used to store one image.
    eHImageFormat ImageFormat;         //See formats below
};
```

```
enum eHCompression
```

```
{
    H_COMPRESSION_NONE=0,
    H_COMPRESSION_JPEG,
    H_COMPRESSION_RLE,
    H_COMPRESSION_HUFFMAN,
    H_COMPRESSION_LZ,
    H_COMPRESSION_RLE_FAST,
    H_COMPRESSION_HUFFMAN_FAST,
    H_COMPRESSION_LZ_FAST,
    H_COMPRESSION_H264,
    H_COMPRESSION_WAVELET
};
```

```
enum eHMetadataFormat
```

```
{
    H_METADATA_UNKNOWN = 0,
    H_METADATA_BOOL,      //bool (1 byte)
    H_METADATA_BYTE,      //byte (1 byte)
    H_METADATA_SHORT,     //short (2 bytes)
    H_METADATA_USHORT,    //unsigned short (2 bytes)
    H_METADATA_INT,       //int (4 bytes)
    H_METADATA_UINT,      //unsigned int (4 bytes)
    H_METADATA_DOUBLE,    //double (8 bytes)
    H_METADATA_STRING,    //wchar_t[] (variable)
    H_METADATA_BINARY     //BYTE[] (variable)
    H_METADATA_INT64,     //__int64 (8 bytes)
    H_METADATA_UINT64     //unsigned __int64 (8 bytes)
};
```


```
enum eHImageFormat
```

```
{
    H_IMAGE_UNKNOWN = 0,           //Unknown format
    H_IMAGE_MONO = 100,            //Monochrome Image (LSB)
    H_IMAGE_MONO_BAYER = 101,     //Raw Bayer Image (treated as H_IMAGE_MONO)
    H_IMAGE_BGR = 200,            //BGR Color Image
    H_IMAGE_PLANAR = 300,         //Planar Color Image
    H_IMAGE_RGB = 400,            //RGB Color Image
    H_IMAGE_BGRx = 500,           //BGRx Color Image
    H_IMAGE_YUV422 = 600,         //YUV422
    H_IMAGE_YUV422_20 = 610,
    H_IMAGE_UVY422 = 700,         //UVY422
    H_IMAGE_UVY411 = 800,         //UVY411
    H_IMAGE_UVY444 = 900,         //UVY444
    H_IMAGE_BGR555_PACKED = 905,  //PhynxRGB
    H_IMAGE_BGR565_PACKED = 906,
    H_IMAGE_MONO_MSB = 112,        //Only for > 8 bit per pixel, MSB align little endian 10 bit: JIHGFEDC BA000000
    H_IMAGE_MONO_BAYER_MSB = 113, //Only for > 8 bit per pixel, MSB align
    H_IMAGE_MONO_MSB_SWAP = 114,  //Only for > 8 bit per pixel, MSB align big endian 10 bit: BA000000 JIHGFEDC
    H_IMAGE_MONO_BAYER_MSB_SWAP = 115, //Only for > 8 bit per pixel, MSB align
    H_IMAGE_BGR10_PPCKED = 123,    //Only for 10 bit per pixel, LSB align
    H_IMAGE_BGR10_PPCKED_PHOENIX = 124, //Only for 10 bit per pixel, LSB align, RRRRRRRR RR00GGGG GGGGGGBB BBBBBBBB
    H_IMAGE_RGB10_PPCKED_PHOENIX = 125, //Only for 10 bit per pixel, LSB align, BBBBBBBB BB00GGGG GGGGGGRR RRRRRRRR
    H_IMAGE_MONO_PPCKED = 131,     //Only for > 8 bit per pixel, MSB align
    H_IMAGE_MONO_BAYER_PPCKED = 132, //Only for > 8 bit per pixel, MSB align
    H_IMAGE_MONO_PPCKED_8448 = 133, //Only for > 8 bit per pixel, MSB align
    H_IMAGE_MONO_BAYER_PPCKED_8448 = 134, //Only for > 8 bit per pixel, MSB align
    H_IMAGE_GVSP_BGR10V1_PACKED = 135, //Only for 10 bit per pixel(From Gige Vision) BBBBBBBB GGGGGGGG RRRRRRRR 00BBGGRR
    H_IMAGE_GVSP_BGR10V2_PACKED = 136, //Only for 10 bit per pixel(From Gige Vision)00BBBBBB BBGGGGGG GGGGRRRR RRRRRRRR
    H_IMAGE_BASLER_VENDOR_SPECIFIC = 1000,
    H_IMAGE_EURESYS_JPEG = 1001,
    H_IMAGE_ISG_JPEG = 1002
};
```


Additionnal Structures and Enums (2/2)

```
enum eHCompression
{
    H_COMPRESSION_NONE=0,           //Uncompressed sequence
    H_COMPRESSION_JPEG=1,           //JPEG compressed images
    H_COMPRESSION_RLE=2,            //RLE compressed images (stp3 algo)
    H_COMPRESSION_HUFFMAN=3,        //HUFFMAN compressed images (stp3 algo)
    H_COMPRESSION_LZ=4,             //LZ compressed images (stp3 algo)
    H_COMPRESSION_RLE_FAST=5,       //RLE Fast compressed images (ippwrapper)
    H_COMPRESSION_HUFFMAN_FAST=6,   //HUFFMAN Fast compressed images (ippwrapper)
    H_COMPRESSION_LZ_FAST=7,        //LZ Fast compressed images (ippwrapper)
    H_COMPRESSION_H264=8,          //H264 compressed images
    H_COMPRESSION_WAVELET=9,        //Wavelet compressed images
    H_COMPRESSION_CREATEC_FLAT_1D = 10, //Flat model predictor using 1 sample for Createc Bayer compressed images
    H_COMPRESSION_CREATEC_FLAT_2D = 11, //Flat model predictor using 2 samples for Createc Bayer compressed images
    H_COMPRESSION_CREATEC_FLAT_BLUE_VECTOR_1D = 12, //Gradient predictor using the blue channel for Createc Bayer compressed images
    H_COMPRESSION_CREATEC_GREEN_AVERAGE = 13, //Gradient predictor using the green channel average for Createc Bayer compressed images
    H_COMPRESSION_CREATEC_2D_LOSSY = 14, //Flat_2D but with a reduced symbol set, Lossy codec with a maximum error of 1 for Createc Bayer compressed images
    H_COMPRESSION_CREATEC_FLAT_2D_SSE = 15 //SSE2 implementation of Flat_2D for Createc Bayer compressed images
};
```


Metadata File Sections

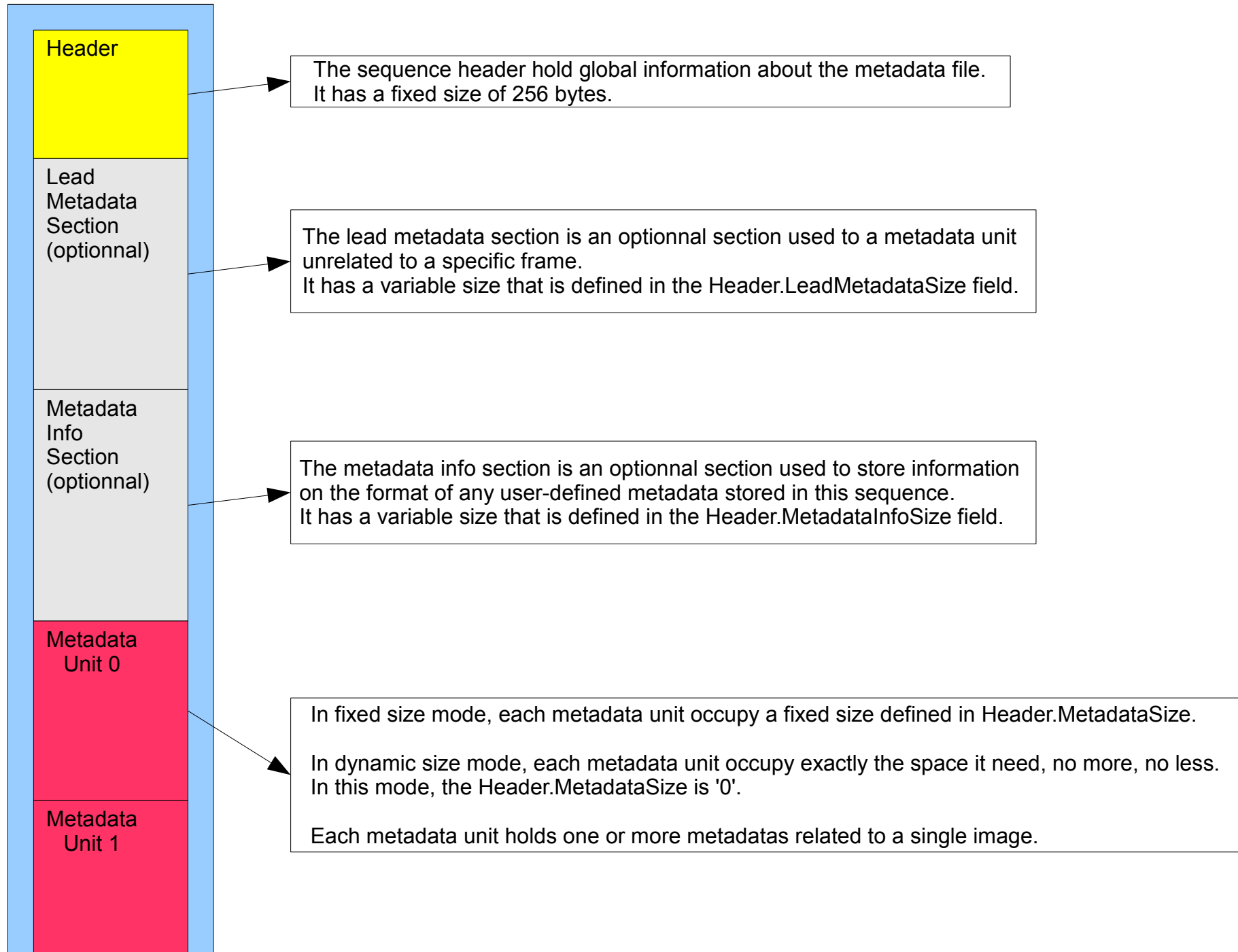


Metadata File Sections

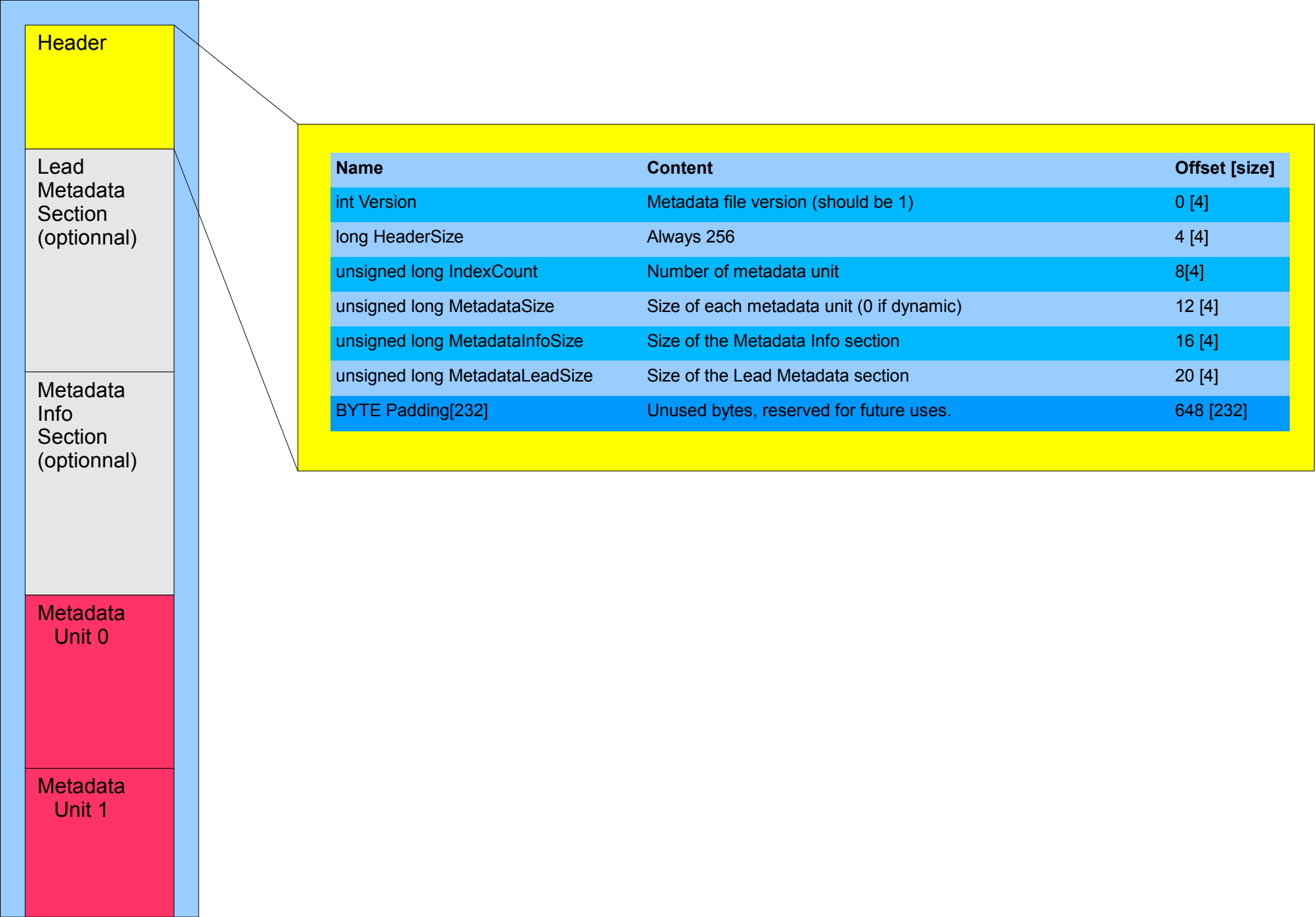

Metadata File Sections

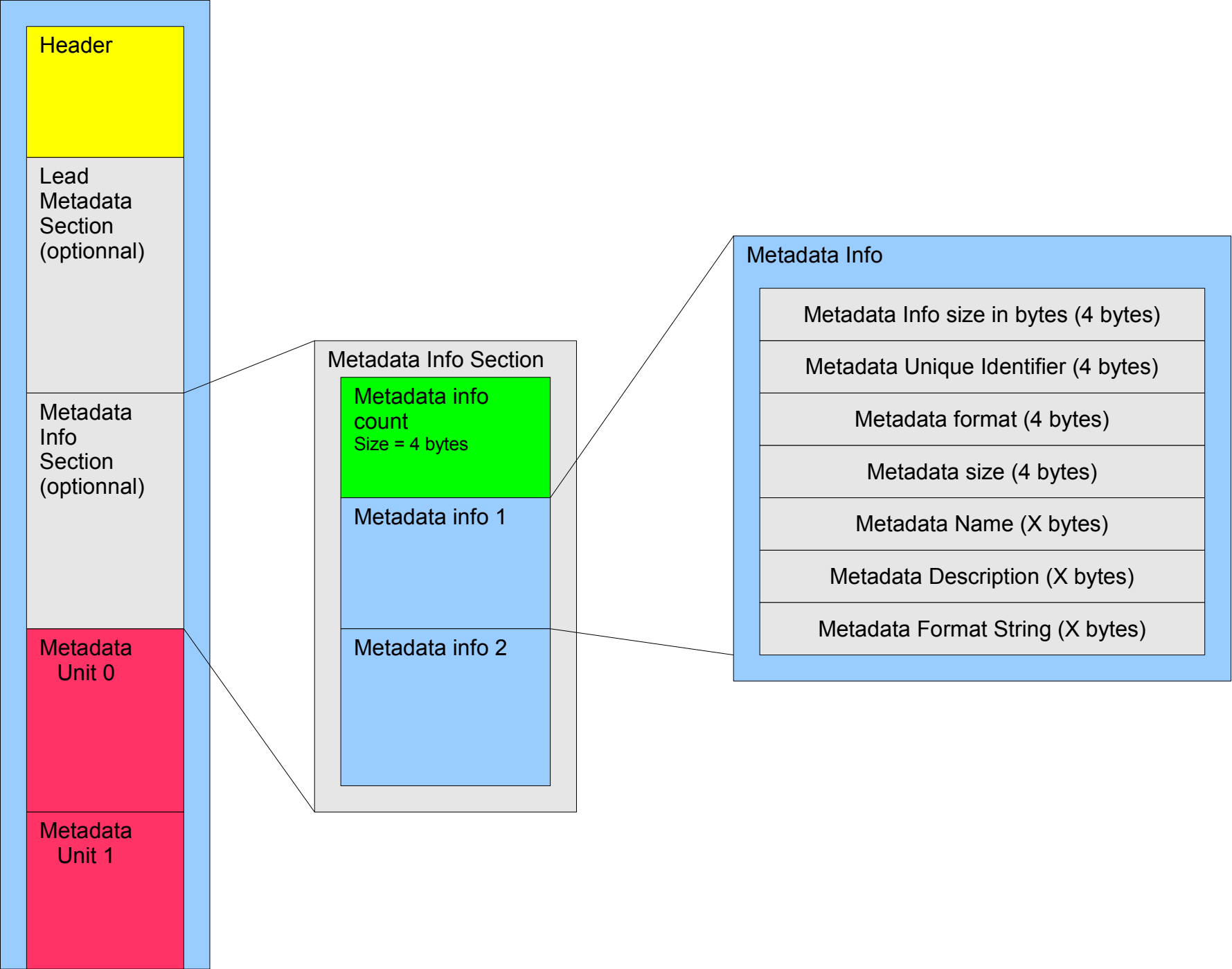

Metadata File Sections

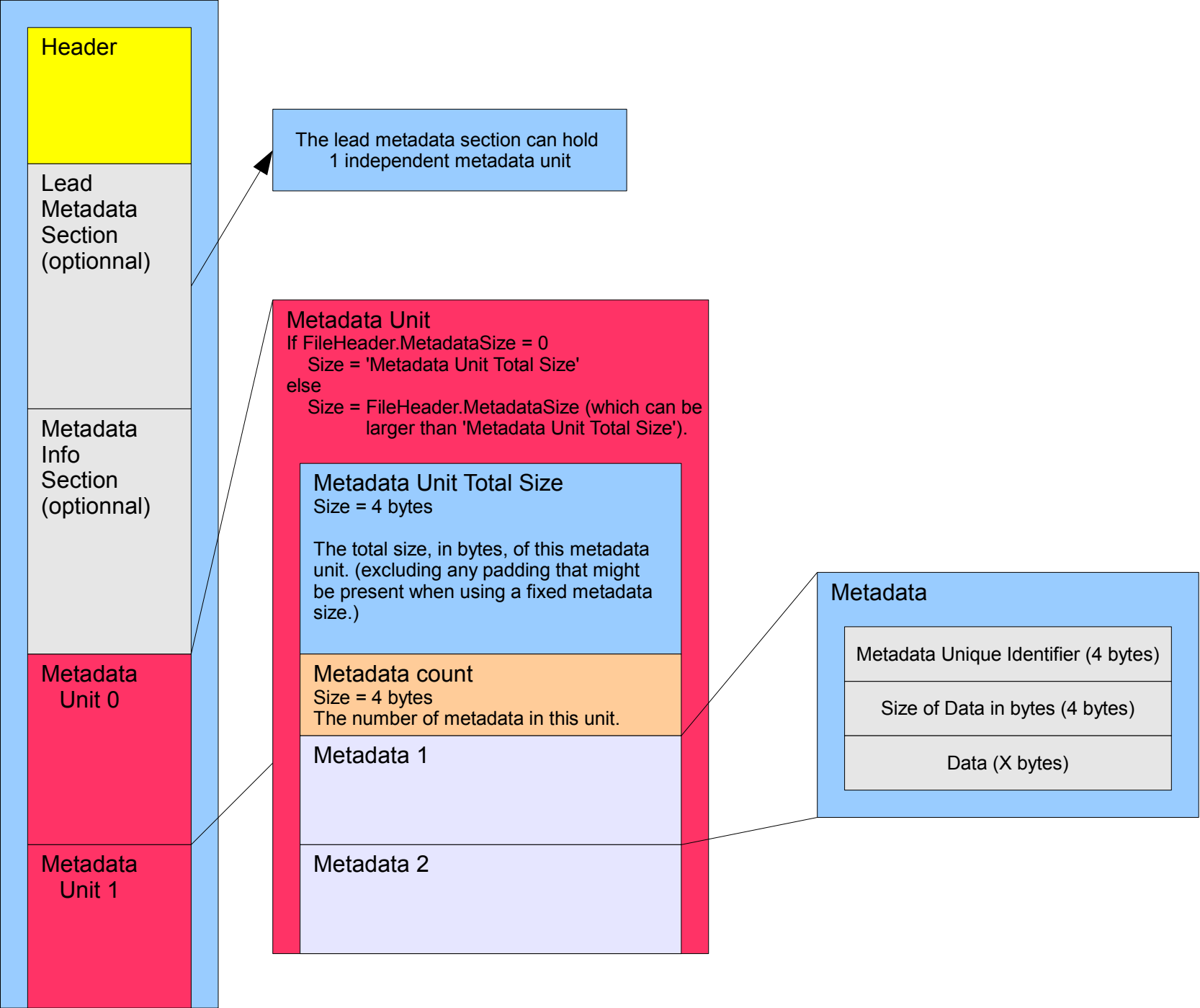

How to set the Metadata size

Fixed size mode

In this mode, the metadata file allocate a fixed size for each metadata unit and a fixed size for the metadata info section. Having a fixed size means that an image's metadata might not be saved to the file if the fixed size is not large enough. It can also lead to waste of space if the fixed size is much larger than what is needed to store the metadata after each image. This mode should only be used when recording in loops.

Dynamic size mode

In the dynamic mode, the metadata file allocate exactly the requiered space needed to store each metadata unit. The metadata info space is allocated by looking at the metadata bundled in the first frame to be written to the sequence. Dynamic mode ensure that you will have exactly the required size to store the metadata (and the metadata info, if any).

Metadata Manager

The Metadata Manager is a stand-alone application that can be used to configure which mode to use and the sizes of the sections. The Metadata Manager also offer a Size Calculator tool that allows you to estimate the required size for the "Metadata" and the "Metadata Info" sections. Simply select the metadata types that are to be saved with each image and the will do the rest.

